


## ATTACHMENTS

Two patented\* attachments are included with the Veinlite EMS PRO:


### Veinlite LED+ Light Shield

A snap-on light shield to reduce interference from bright overhead lighting.

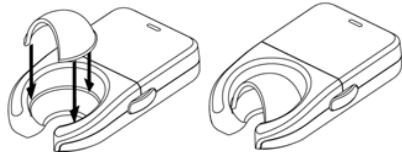

### Veinlite LED+ Pediatric Adapter

A snap-on pediatric adapter to reduce the size of the viewing area for use with infants.

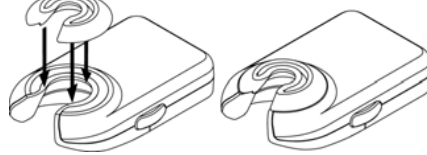




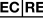

## WARNINGS

- Veinlite LED+ is designed for external examination only.
- To avoid contamination of the device, Veinlite LED+ must be used with the disposable plastic cover.
- To avoid the spread of infection, always change the disposable plastic cover between patients and wipe the new cover with a germicidal wipe.
- These devices are intended for use only by clinicians trained in performing vascular access, or under the supervision of a clinician.
- Do not look directly into the LEDs while on.

## Warranty Information

The Veinlite LED+ has a one year parts and labor warranty. Please provide the serial number located under the battery or the outside of the box when you contact TransLite LLC for any warranty issues. The warranty is valid only if the Veinlite LED+ is used according to the manufacturer's instructions. Before returning the Veinlite LED+ to TransLite for repair, you must call or e-mail TransLite to obtain a return authorization code and shipping information before sending your device. No refunds will be given for unauthorized shipping.

## Explanation of Symbols

-  Model Number
-  Serial Number
-  Manufacturer
-  EC Representative (EU Only)
-  Do not dispose of as unsorted municipal waste


 **TransLite LLC**  
345 Commerce Green Blvd  
Sugar Land TX 77478 USA  
+1 281 240 3111  
www.veinlite.com  
info@veinlite.com

 **M. Devices Group**  
Health Education Centre  
The Church, Portland Street  
Southport PR8 1HU  
United Kingdom  
+44 1704 544 944

\* This product covered by US Patents:  
7,874,698; 7,841,751; 8,032,205; 8,177,808

© 2018 TransLite, LLC. All rights reserved.  
These instructions or any part thereof may not be reproduced by any means, without prior permission of TransLite, LLC.

**Made in USA**   
VLEDP DFU Rev 4 0220

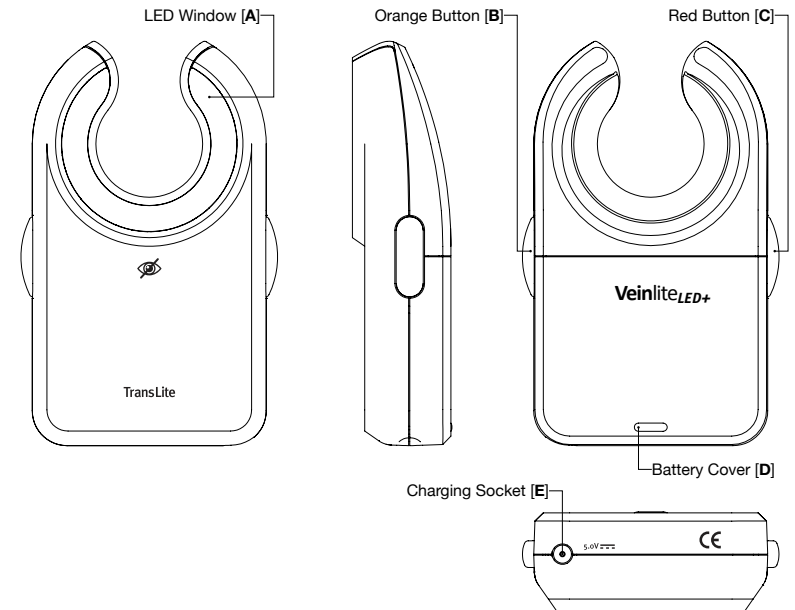

## DIRECTIONS FOR USE

### Description

The Veinlite® LED+ is a hand-held battery powered transillumination device for visualizing veins and superficial blood vessels. Veinlite utilizes the patented\* transillumination technique known as side-transillumination, which uniformly illuminates a small region of skin and subcutaneous tissue without any shadow areas. The shadow-free side-transillumination technique allows better visualization of veins than any other transillumination method.

The LEDs emit light with the optimal wavelength for highlighting the contrast between the venous blood vessels and surrounding tissue.

### Intended Use/Indications

- Venous access in children/adult patients
- Patients with a history of difficult venous access
- Improved visualization of veins in patients with different pigmented skin types
- Improved visualization of veins in geriatric patients

### Contraindications


- When used for internal examination
- When used without the disposable plastic cover




## OPERATION

- Place the Veinlite LED+ against the patient's skin with the LED Window [A] in direct contact with the skin perpendicular to vein direction. Apply gentle pressure to maintain optimal skin contact.
- To turn on the LEDs, press the Orange Button [B] or Red Button [C] located on the left and right side of the device.
- Press the button a second time to turn off.
- Press the Orange Button [B] if you are examining small superficial veins on fair skin.
- Press the Red Button [C] if you are examining deeper veins or if the patient's skin is darkly pigmented.


**WARNING: DO NOT SHINE THE LED LIGHTS IN THE PATIENT'S EYES.**

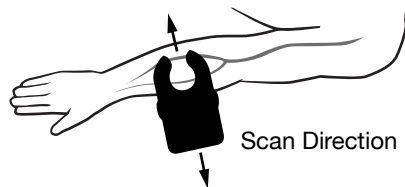


- 1 To scan for veins, place device against skin perpendicular to vein direction. Scan it forward and backward. Veins will be visible as dark lines due to transillumination.

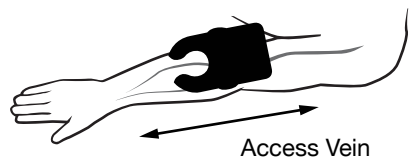

- 2 Once a suitable vein has been identified, rotate device 90°, parallel to vein.

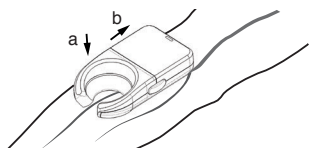

- 3a Apply gentle pressure on the Veinlite LED+ to secure the vein.
- 3b Push back slightly on the Veinlite LED+ to stretch the skin and provide traction.

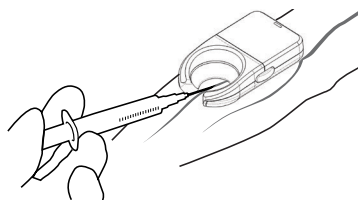

- 4 Insert the needle into the vein, through the opening of the device.  
NOTE: Discard the single-use disposable plastic cover after each patient.

## USING THE DISPOSABLE PLASTIC COVER

The single-use disposable plastic covers for the Veinlite LED+ (Product #VLED+DPC) MUST be used when the device is placed in contact with a patient's skin.

CAUTION: These plastic covers are not sterile. They should be wiped with a germicidal wipe prior to use.

- Before using Veinlite LED+, place it inside the disposable cover with the flap on the same side as the LEDs.
- Pull the flap over the device to the other side.
- Wipe the LED side of the Veinlite LED+ and any part of it coming in contact with the patient, with a medically approved germicidal wipe.

## PRECAUTIONS

- Only use device with the accompanied charging cord supplied by TransLite.
- Do not use the device while plugged into any power source.
- Device may be wiped down with a germicidal wipe however, avoid using alcohol on the clear plastic window as this may cause fogging.
- Do not immerse the device in liquid.
- Do not autoclave the device.

### Battery Level Indicator Light

Veinlite LED+ incorporates a low battery level indicator light. The orange button illuminates when the battery charge level has fallen to approximately 5% of the full battery charge.

### Recharging the Battery

The Veinlite LED+ is powered by a rechargeable lithium ion battery, which is designed to provide approximately 300 minutes of continuous usage.

- Connect the charger to the Veinlite LED+ by inserting it into the Charging Socket [E] located on the bottom of the device.
- The orange button will blink indicating that the device is charging.
- Once the blinking stops your device is fully charged.

### Replacing the Battery

- Slide off the Battery Cover [D].
- Remove the battery by unplugging the connector at the top.
- Replace the battery with a new battery by plugging it into the connector
- Slide the Battery Cover [D] back on.

NOTE: Before using your Veinlite LED+, please read the detailed battery safety information on our website FAQ page: <https://www.veinlite.com/faq>.